

 rouge

**User
Manual**




What is Red Therapy Light and How Does it Work?

Red light therapy is also known as low-level laser therapy (LLLT). It uses low-energy light emitting diodes (LEDs) or lasers as an alternative therapy for pain relief, or to promote good cell functions. The principle of red light therapy can be thought of as like a human version of the way plants photosynthesize. It is thought that certain

wavelengths of light can enhance the functioning of mitochondria, which are the organelles in a cell which take in nutrients and break them down to create energy molecules for the cell. By enhancing the cellular respiration, it will improve the functioning of the whole cell.


Benefits Of Rouge LED Light Therapy

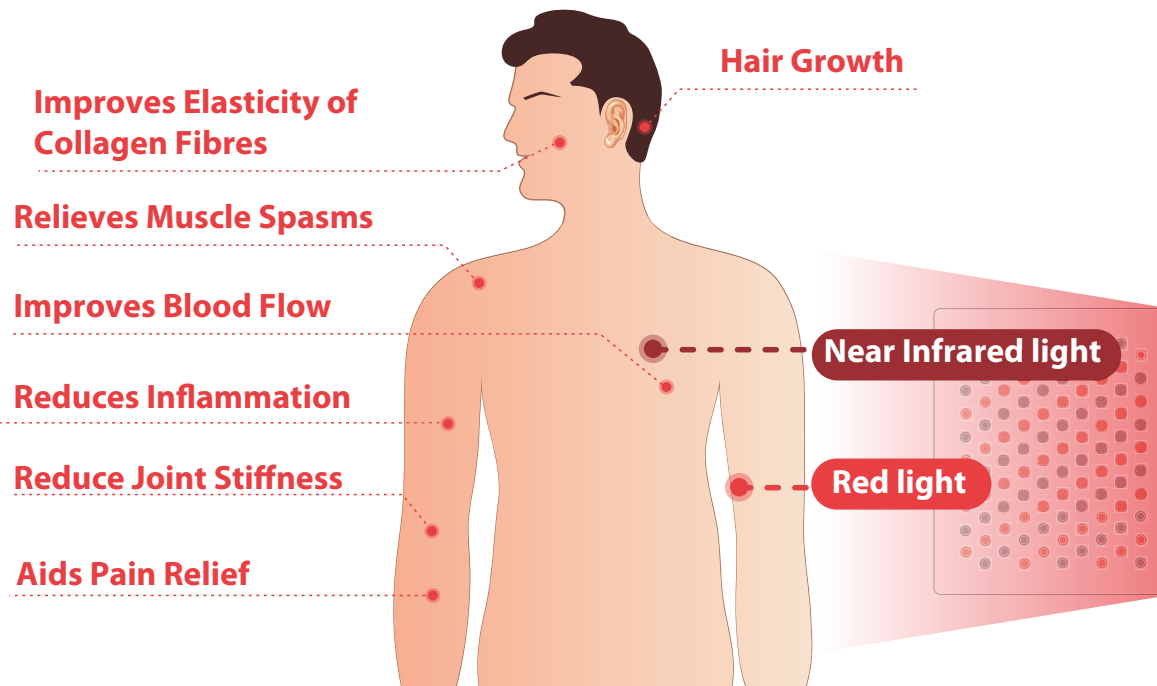



Why choose us?


1. **More Power Density:** Independent tests show Rouge is in a class of its own. With 2.75X more power density than our main rival, Joovv, it's no surprise that Rouge has so quickly become the #1 choice of naturopaths and serious biohackers, worldwide.
2. FDA registration is just the beginning - our products have received numerous safety certifications and are among the world's safest and most reliable red light therapy products.
3. **Satisfaction Guaranteed:** If for any reason you're not totally satisfied with your panel, you can return it to us within 60 days for a full refund.

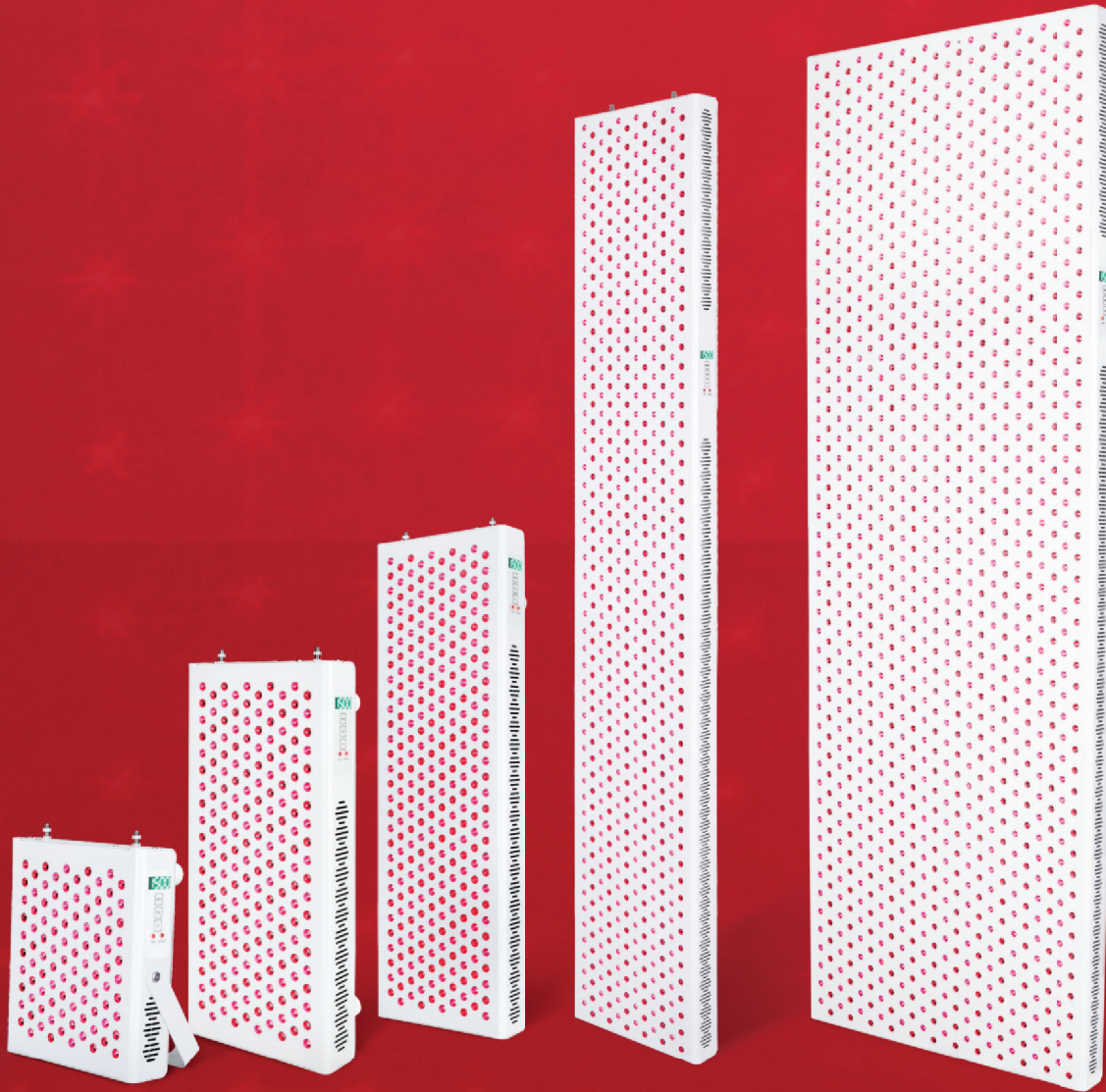


Rouge Product Specifications

	X	TABLETOP	ESSENTIAL	PRO	MAX	ULTIMATE
Wavelengths:	RED 660nm & NIR 850nm	RED 660nm & NIR 850nm	RED 660nm & NIR 850nm	RED 660nm & NIR 850nm	RED 660nm & NIR 850nm	RED 660nm & NIR 850nm
Dimensions	18.1" x 13.8" x 50.4"	10.4" x 8.7" x 2.3"	18.9" x 8.7" x 2.3"	27.2" x 8.7" x 2.3"	64.9" x 11.8" x 2.3"	65" x 21.7" x 2.3"
Weight	5KG/11LBS	2.8KG/6.2LBS	4.5KG/9.9LBS	7KG/15.4LBS	21KG/46.3LBS	36KG/79.4LBS
Irradiance (mW/cm2)*	130	135	145	161	103	100
EMF Emission	- Negligible	Negligible	- Negligible	- Negligible	- Negligible	- Negligible
Amps @110v	1.15A	1.1A	2.1A	3.3A	8.68A	13.4A
Amps @220v	0.68A	0.4A	1.1A	1.3A	2.92A	6.83A
Power Consumption	82W	120W	235W	352W	910W	1443W
Beam Angle	30 Degree	30 Degree	30 Degree	30 Degree	30 Degree	30 Degree
Best For	Targeted Treatment	Targeted Treatment	Targeted Treatment	Full Body Treatment	Full Body Treatment	Full Body Treatment
Door Mount	✓	✓	✓	✓	✗	✓
Built-in Programmable Timer	✓	✓	✓	✓	✓	✓
Smart Built-in Cooling Fans	✓	✓	✓	✓	✓	✓


Irradiance Charts

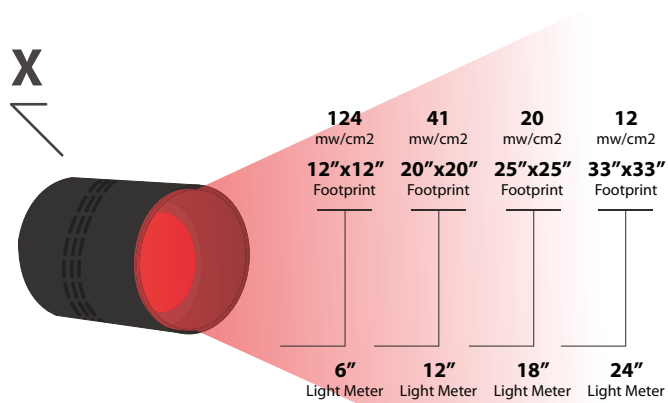

ESSENTIAL

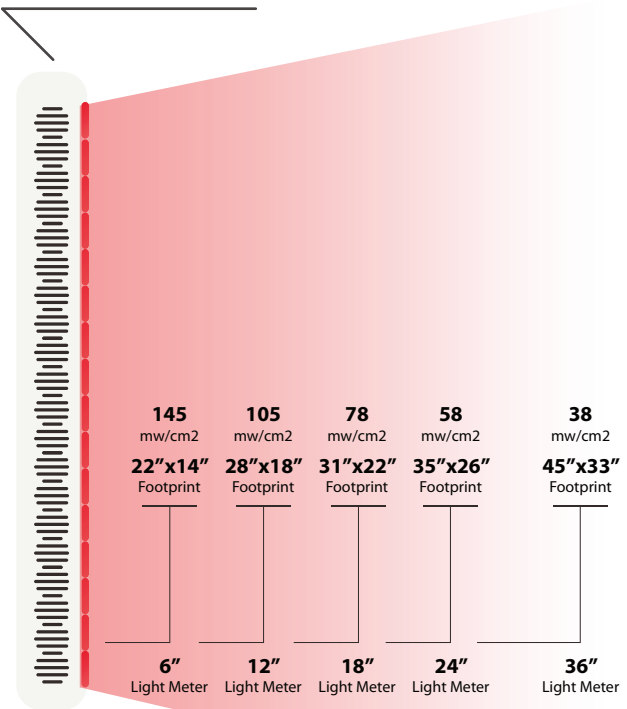


TABLETOP

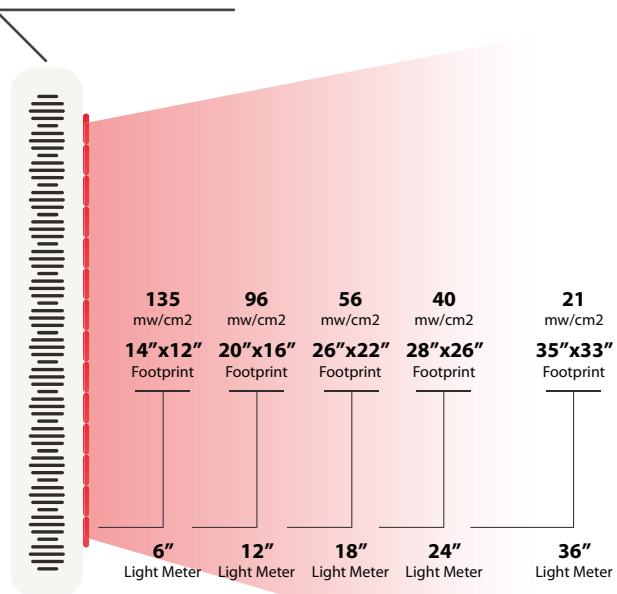


Irradiance Charts


PRO

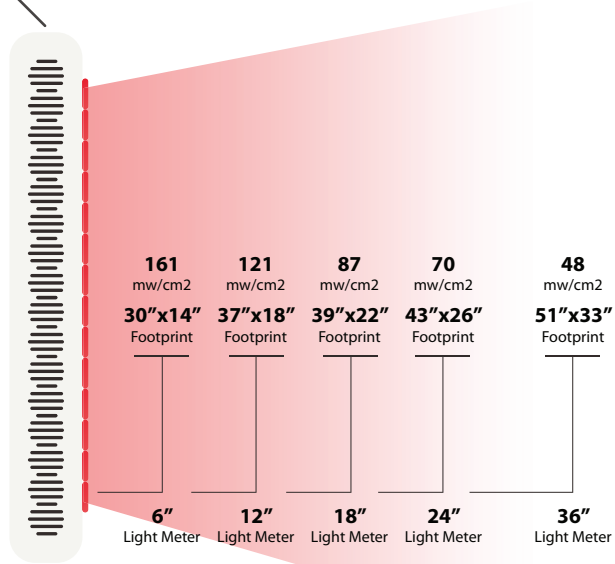

MAX

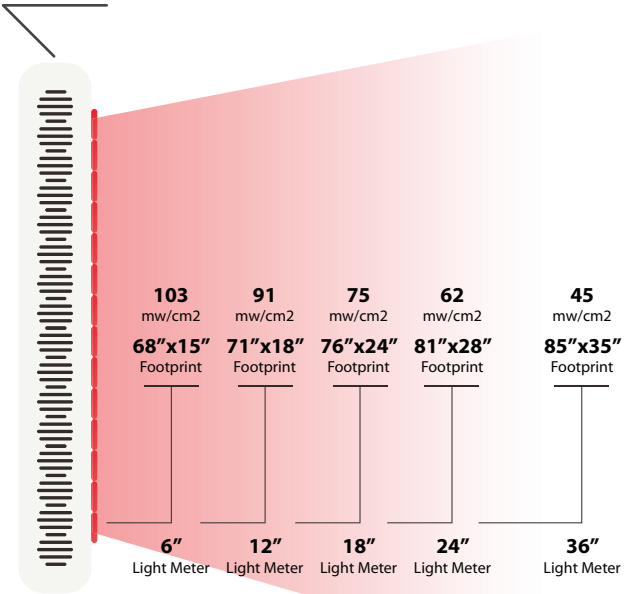

ULTIMATE

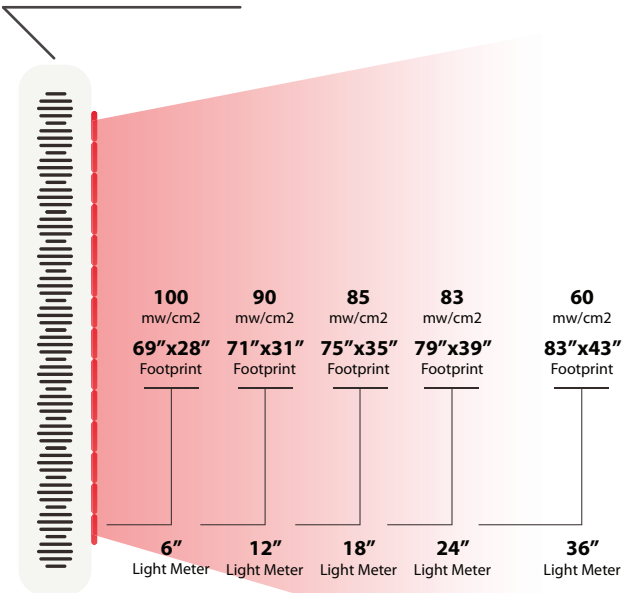


Optional Mounting Configurations



Tabletop Stand



Door Mount

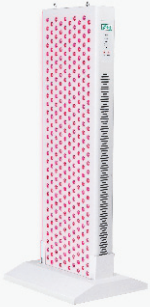

Floor Stand



Wall Setup



Wheeled Standing Rack



Package Contents

Included with all
Devices:

- 1 Rouge LED Device
- 4 Power Cable
- 5 Goggles

Included only with
Tabletop, Essential & Pro:

- 2 Patch Cable
- 3 Pulley Assembly & Snap Link
- 6 Hanging Cables
- 7 Door Hook


Included only with
Ultimate & MAX:

- 8 Wall Mount Bracket

Install Your Therapy Light


STEP 1 ▶


Lay all of the included hardware, including your Rouge LED, on a padded surface. Next, screw in the threaded posts to the top of the device and hand-tighten with the included Allen wrench. Then, screw the cable supports to the threaded posts, as shown in PIC.1.

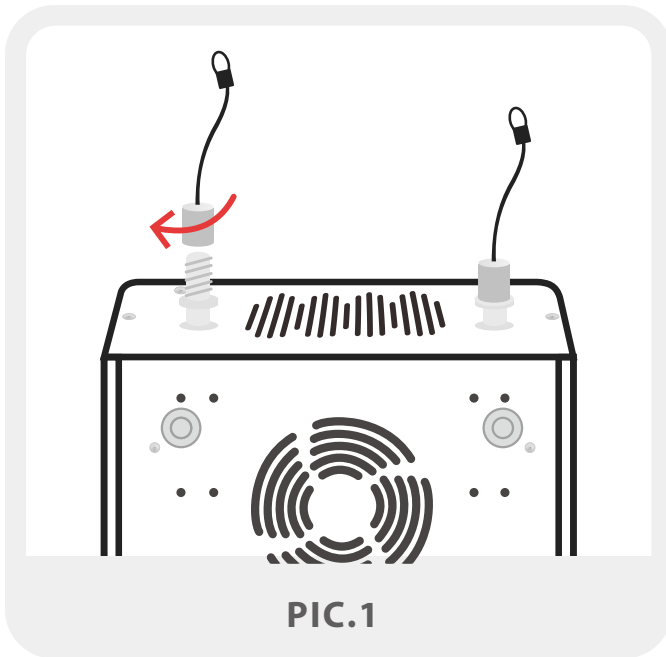

STEP 2 ▶


Hang the included custom hook on any standard-sized door, as shown in PIC.2.

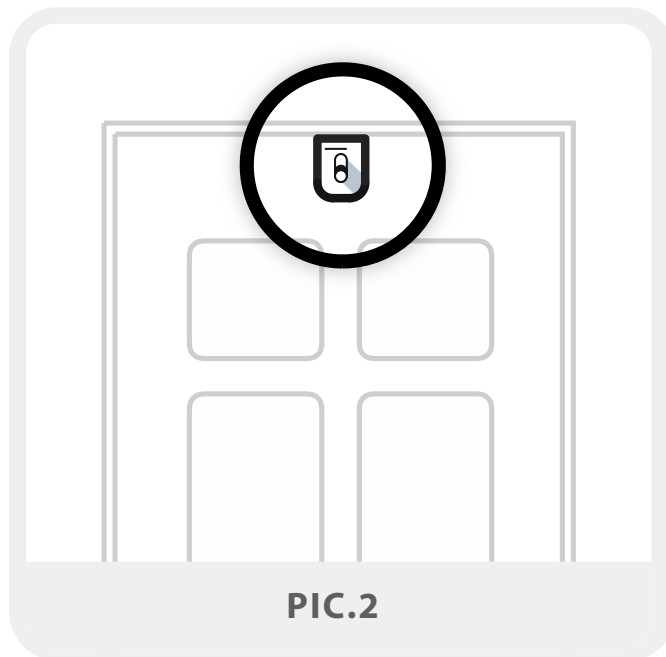

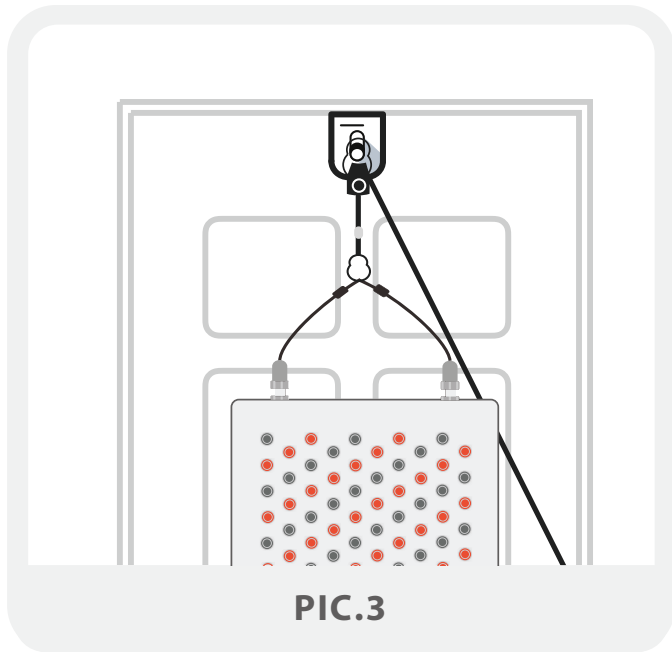

STEP 3 ▶

Connect the cable supports to the snap link attached to the height-adjusting system. Then, connect the snap link on the opposite end directly to the door hook, as shown in PIC.3.

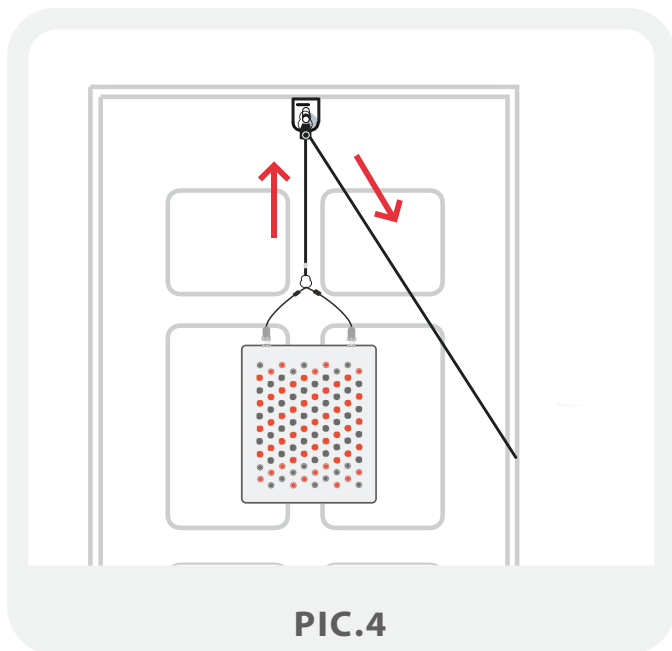


STEP 4 ▶

To raise your RougeLED to your desired height, pull down on the loose cable while supporting the device with your opposite hand, as shown in PIC.4.




Be sure to adequately support the device while raising and lowering it with the height-adjusting system .

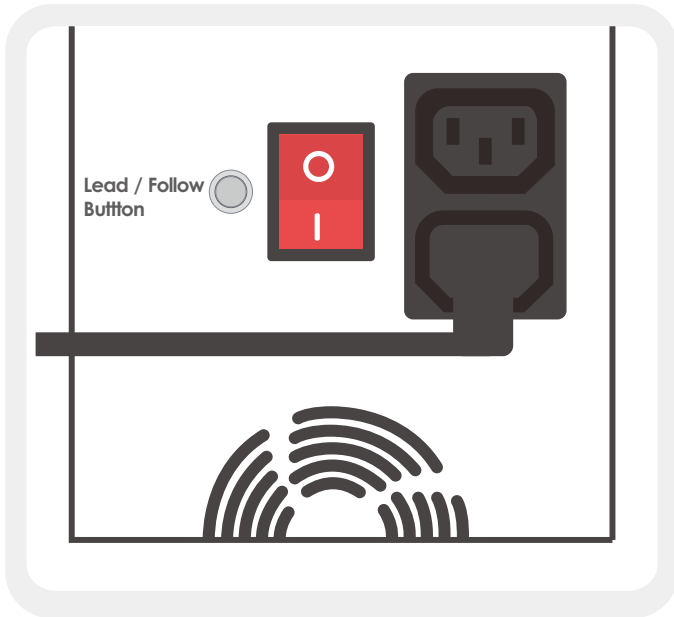

Setting up Your Rouge Pro Series Light

STEP 1 ▶

Plug the power cord in the female socket and press the red switch, it will become illuminated.

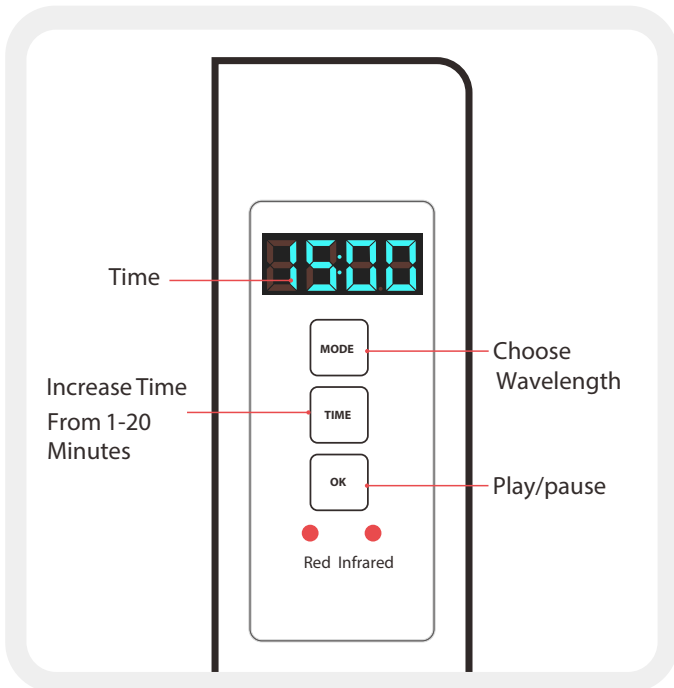


STEP 2 ▶

1. You can choose only red or only IR or both by selecting the "Mode" button. The red indicator light will illuminate indicating which wavelength configuration is selected.
2. Adjust the time by pressing the "Time" button.
3. Press "ok" to run the light, start the treatment.




Remember that near infrared LEDs emit light that is NOT visible to the naked eye.

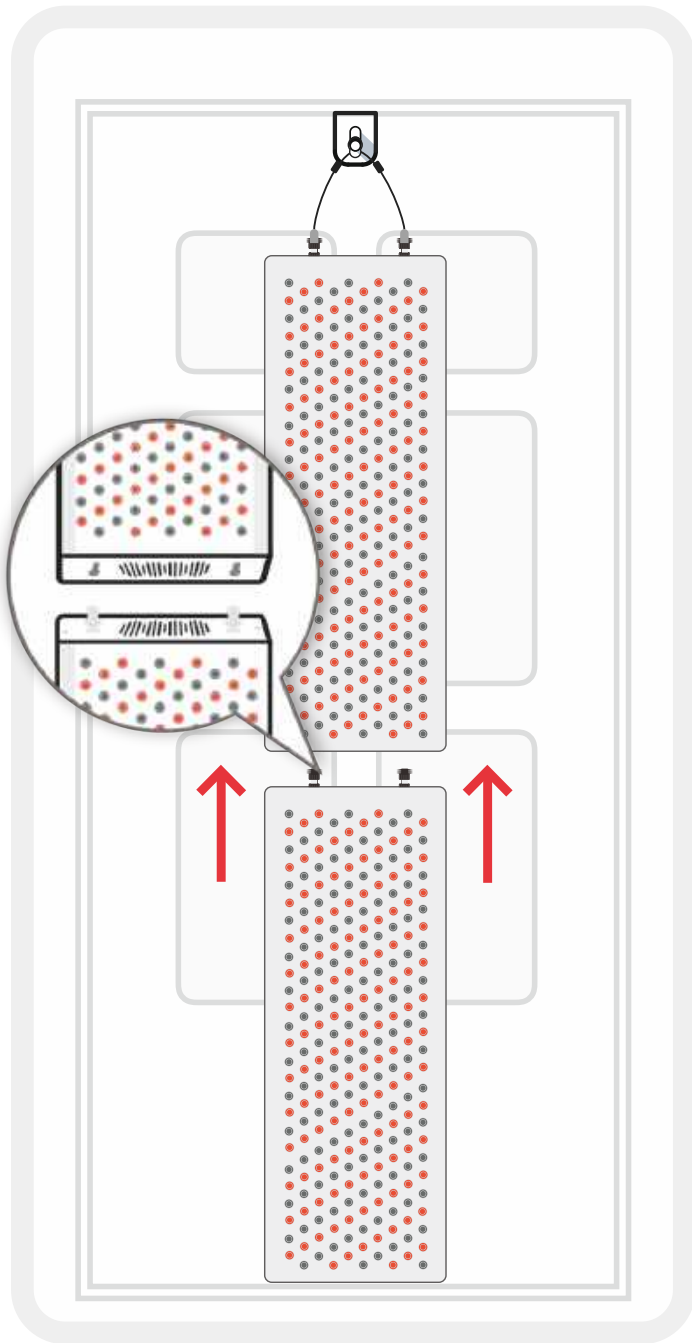

Connecting Multiple Rouge Pro Series Light Units

STEP 1 ▶

To connect the lower Rouge LED, lift it into position below the lead device that is already in place. Line up the threaded posts with the holes in the end caps of the upper device. Then, slide the lower threaded posts into the holes, and slide forward to lock in place.

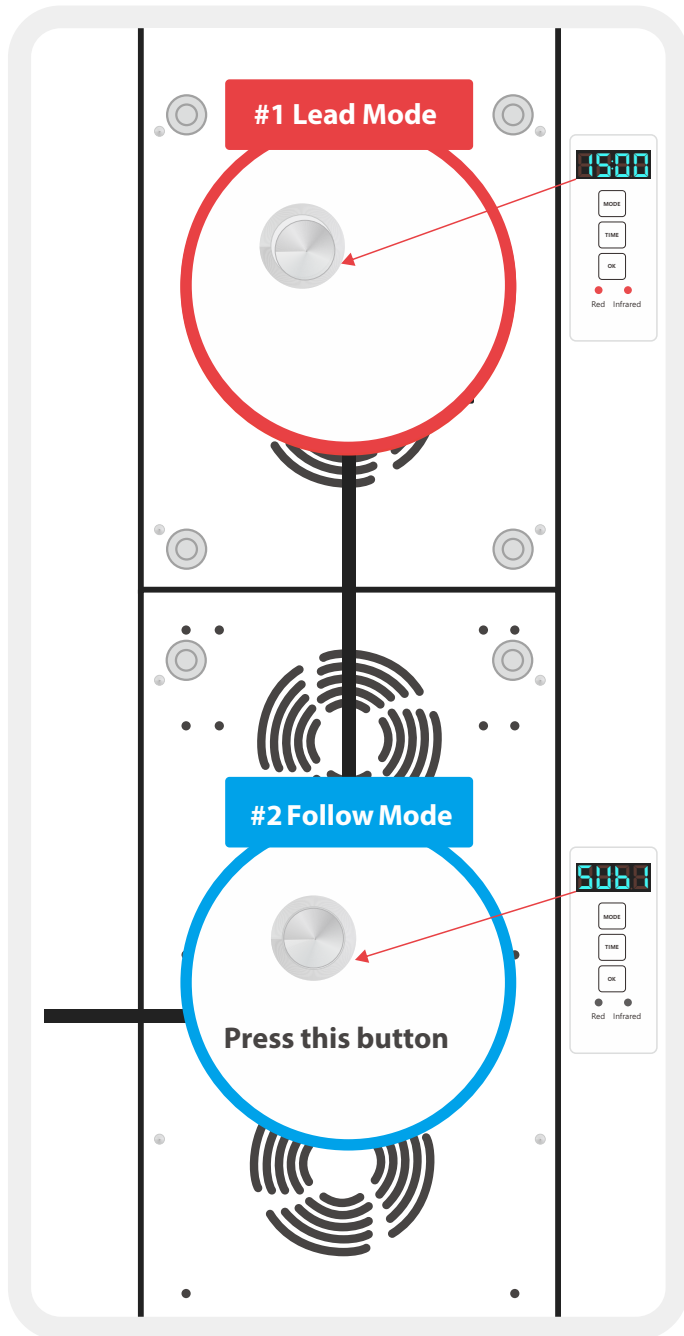

STEP 2 ▶

1. Choose #1 light as leader.
2. Press #2 light's follow button, the screen will show "SUB1", #2 lights stay in "Follow Mode".
3. You only need to control the lead light, other lights will follow the lead light's instruction.

Rouge Warranty

Warranty Clauses

1. Rouge Red Light Therapy devices are warranted against defects in materials and workmanship for a period of THREE (3) YEARS from the purchase date.
2. Free replacements are warranted for a period of SIX (6) MONTHS.
3. Warranty does not cover the following damage(s): Modified product, improper installation or use, improper cleaning, electrical surges, or damage caused by direct exposure to natural elements like water or fire for example.
4. Warranty does not cover the loss of third party with direct or indirect damage(s)

Warning & Cleaning

1. Unplug fixture before cleaning.
2. Do not use cleaning sprays. Wipe with a dry or lightly damp cloth to clean.



rouge.care

✉ E-Mail: info@rouge.care

☎ Tel: 1 (833) 848-2977

👤 Rouge Red Light Therapy
6000 Boul Rome, suite 410 Brossard, QC J4Y 0B6

©2020 -2021 Rouge Red Light Therapy. All rights Reserved